



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
STATE PUBLIC HEALTH LABORATORY
BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

RECEIVED

By Carol Day at 8:39 am, Oct 21, 2013 #6

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).
Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.
Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 201083	NAME OF AGENCY Boone Co. S.O.	DATE OF INSPECTION 10/20/2013
LOCATION OF INSTRUMENT (STREET AND CITY) 2121 County Dr. Columbia		TIME OF INSPECTION 19:15

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 10/20/13 19:15
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 49 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER

<input checked="" type="checkbox"/> INDICATOR LIGHTS
<input checked="" type="checkbox"/> SIMULATOR SOLUTION SUPPLIER Guth Labs LOT # 13100 EXP. DATE 04/23/2015
<input checked="" type="checkbox"/> SIMULATOR TEMP (34°C ± 0.2°C) 34.0 °C SIMULATOR SN SD2218 EXP. DATE 06/14/2014

<input checked="" type="checkbox"/> CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)		
Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)		
<input checked="" type="checkbox"/> 0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE		
<input type="checkbox"/> 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE		
<input type="checkbox"/> 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE		

TEST 1 •100	TEST 2 •100	TEST 3 •101
--------------------	--------------------	--------------------

<input checked="" type="checkbox"/> PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT: (DO NOT INCLUDE SELF-ADMINISTERED TESTS)		
--	--	--

REFUSALS 0	(0-.04) 0	(.05-.09) 2	(.10-.14) 5	(.15-.19) 1	OVER .19 1
-------------------	------------------	--------------------	--------------------	--------------------	-------------------

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER	
SIGNATURE 	PRINT FULL NAME Domenica P. Antini
TYPE II PERMIT NUMBER/EXPIRATION DATE 220294 09/25/2014	TELEPHONE NUMBER 573-875-1111
RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office 2875 James Blvd. Poplar Bluff, MO 63901	

BAC DataMaster
Evidence Ticket

STATE OF MISSOURI
BAC DATAMASTER SERIAL NUMBER 281083
10/20/13

ARREST TIME: 18:00

SUBJECT NAME:

TEST/RFI

I.D.#: 81/81/10 SEX: F

STATE/D.L.: MO/123456

ARRESTING OFFICER:

ANTIMI/N

OFFICER I.D.: 5283

TESTING OFFICER:

ANTIMI/NIKKI

OFFICER I.D.: 5283

PERMIT NUMBER: 220294

EXPIRATION DATE: 09/25/14

MISCELLANEOUS DATA:

RFI TEST

— BREATH ANALYSIS —

RADIO INTERFERENCE



Operator Signature

2208-02

Face This Side Down - This Edge In First

Face This Side Down - This Edge In First

BAC DataMaster Evidence Ticket

STATE OF MISSOURI
BAC DATAMASTER SERIAL NUMBER 2010683
10/29/13

TESTING OFFICER:
ANTHILL/MIKKI
OFFICER I.D.: 5283
PERMIT NUMBER: 220294
EXPIRATION DATE: 09/25/14
MISCELLANEOUS DATA:

— SUPERVISOR MODE —

BLANK TEST .000 19:17
INTERNAL STANDARD VERIFIED 19:17
EXTERNAL STANDARD .100 19:18
BLANK TEST .000 19:19
EXTERNAL STANDARD .100 19:19
BLANK TEST .000 19:20
EXTERNAL STANDARD .101 19:20
BLANK TEST .000 19:21

 N = 3
 SIM. = .1
 PWS. = .10003

STATE OF MISSOURI
BAC DATAMASTER SERIAL NUMBER 2010683
10/26/13
19:15

— DIAGNOSTIC CHECK —

COMPUTER: OKAY
PROGRAM (04-07-2009): OKAY

HEATERS SAMPLE CHAMBER: 49°C

FLOW DETECTOR: OKAY
PUMP HIGH SPEED: OKAY
DETECTOR: OKAY
FILTERS: OKAY

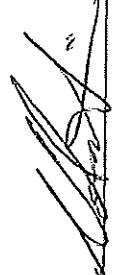
 QUARTZ STABILIZER: OKAY

 CALIBRATION: OKAY

PRINTER TEST
!"#\$/%&'()**,-./0123456789!;=<>?@BCDEFG
HIJKLMNOPQRSTUVWXYZ\^_`abcdefghijklmnopqrstuvwxyz

Operator Signature 

2208-02

Operator Signature 

2208-02

Bottle #601



GUTH LABORATORIES, INC.

892 NORTH 6TH STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13100 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on April 29, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1214% (w/vol) ethyl alcohol. The expiration date for this lot number is April 23, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- 2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:
Testing was conducted using Cerilliant Reference Standard lot number FN1222J1-02 whose values are traceable to NIST.
All balances are calibrated annually by an outside agency using NIST traceable weights.
Calibration verification is done prior to each use utilizing NIST traceable weights.

State of Missouri
DEPARTMENT OF HEALTH

P E R M I T
TYPE II



DOMENICA PANTINI

Is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field repairs, and operate the following breath analyzer(s):

DATAMASTER; ALCO-SENSOR IV W/PRINTER

for the determination of the alcohol content of blood from a sample of expired (alveolar) air. Issued under the provisions of sections 677.020 through 677.041, RSMo 1980.

Date 09/25/2012

Number 220294

Expires 09/25/2014

MO 100-0711 (7-9)

Director of State Public Health Laboratory

Margaret T. Danaher

Director, Department of Health

Lab. # (PZ-K)