

STATE PUBLIC HEALTH LABORATORY
BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

REPORT #8

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days). Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.	
DATAMASTER SN 259 NAME OF AGENCY BOON VILLE PL	DATE OF INSPECTION 08/08/2013
LOCATION OF INSTRUMENT (STREET AND CITY) 401 E. Margan St. BOONVIII	TIME OF INSPECTION O 14
CHECKLIST: Place a mark in the box by each item if found to be sa	alisfactory or if operating within established limits. (Write in observed values
where determined.) Unmarked items must be corrected before us DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) <u>Q8 /8 / 13 2019</u>
	DETECTOR
BPROGRAM (04-07-2009)	₽ FILTERS
HEATERS SAMPLE CHAMBERC	Quartz standard
K FLOW DETECTOR	CALIBRATION
DEPUMP HIGH SPEED	PRINTER
INDICATOR LIGHTS	
2 SIMULATOR SOLUTION SUPPLIER Guth Lab	S LOT # 13/00 EXP, DATE \$ 4/23/2015
SIMULATOR TEMP (34°C \pm 0.2°C) 34.0 °C SIMULATOR SN 50.3501 EXP. DATE $8/8/20/9$	
☑ CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)	
Figure three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of ,005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)	
0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE	
0.040% STANDARD - MUST READ BETWEEN 0.038% A	ND 0.042% INCLUSIVE
TEST 1 096 TEST 2 0	97 TEST 3 = -097
PERFORM R.F.I. TEST (PRINTOUT ATTACHED)	
INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT: (DO NOT INCLUDE SELF-ADMINISTERED TESTS)	
REFUSALS (004) (.0509)	(.1014) / (.1519) OVER .19
LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).	
·	
INSPECTING OFFICER	
SIGNATURE	PRINT FULL NAME E/I'C Lee MOSS
TYPE II PERMIT NUMBER/EXPIRATION DATE	E/I'C Lee 12055
220182 8/03/2019	660 882 2727
RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office 2875 James Blvd. Ponlar Riuff MO 83901	



CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13100 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on April 29, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1214% (w/vol) ethyl alcohol. The expiration date for this lot number is April 23, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.


The alcohol and water used in this solution were free of test interfering substances.


Ted L. Pauley, President GUTH LABORATORIES, INC.


NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

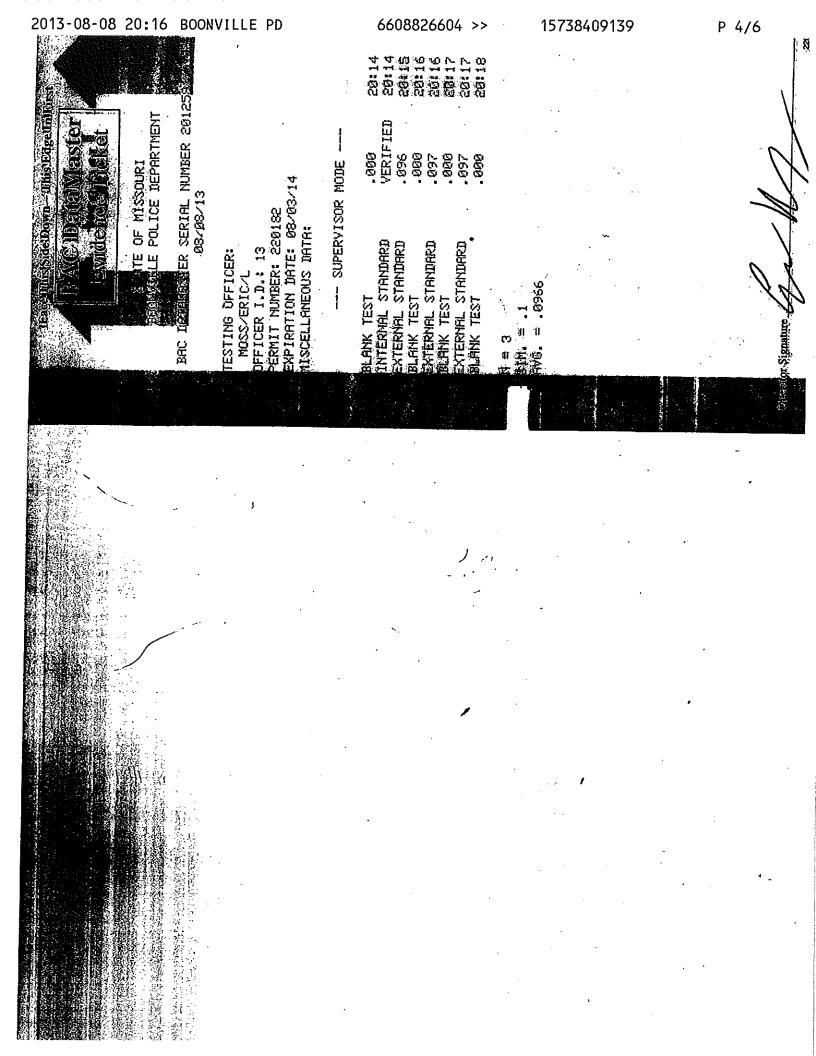

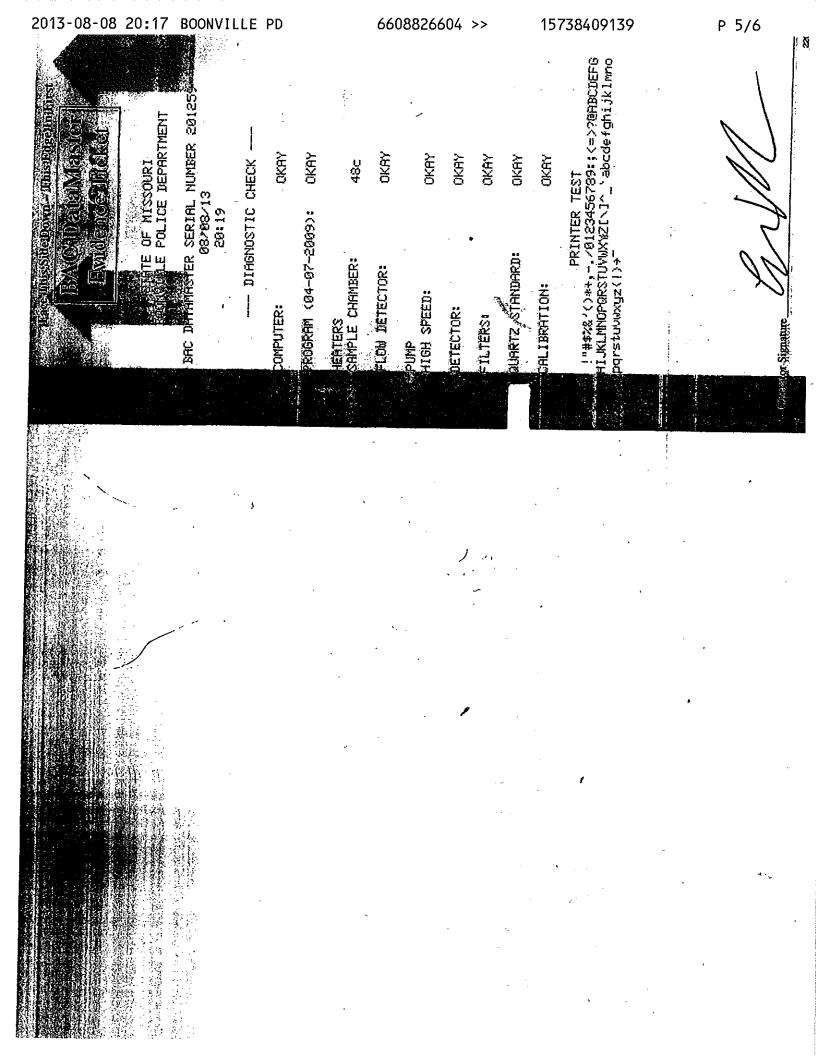

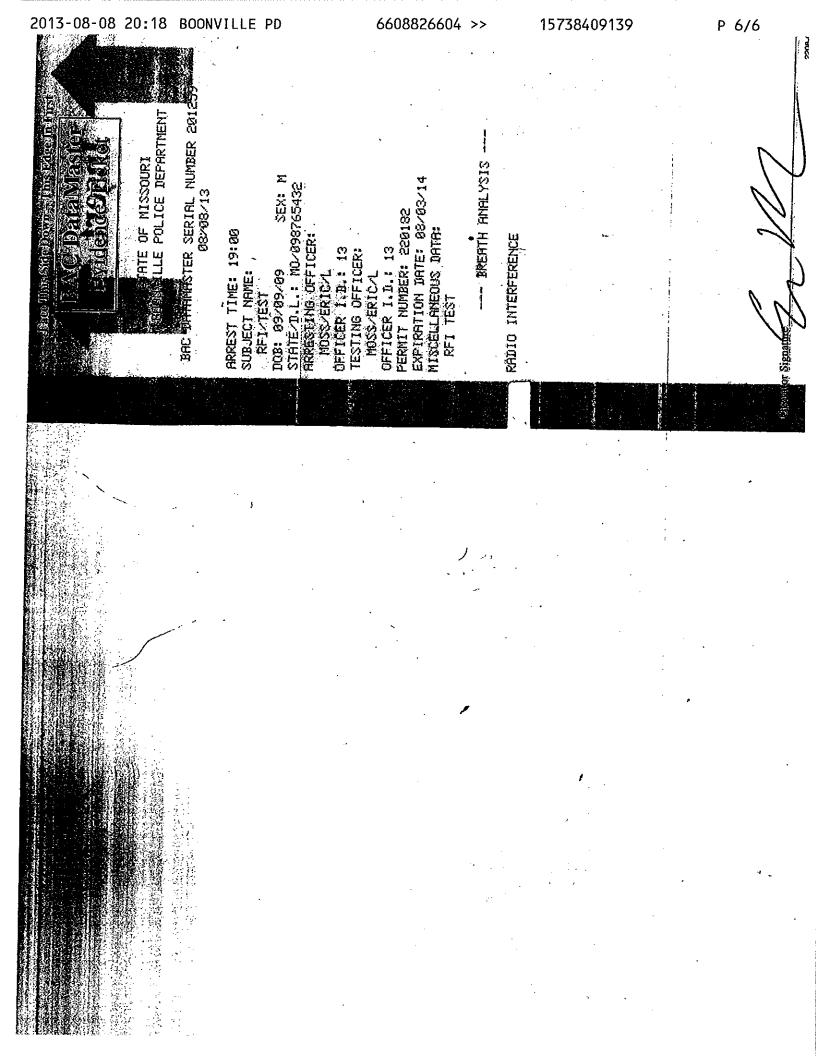

State of Missouri DEPARTMENT OF HEALTH



PERMIT TYPE II



ERIC MOSS

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired (alveolar) air. Issued under the provisions of sections 577.020 through 577.041, RSMo 1986.

08/03/2012

umbar 220182

Expires 08/03/2014

MO 580-0771 (7-88)

Director of State Public Health Laboratory

Director, Department of Health

Lab. 4 (R7-88)