



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES  
STATE PUBLIC HEALTH LABORATORY  
BREATH ALCOHOL PROGRAM

RECEIVED

By Carol Day at 8:23 am, May 28, 2014

INTOX EC/IR II MAINTENANCE REPORT

REPORT #3

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days). Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service. Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

INTOX EC/IR II SN 12710	NAME OF AGENCY St. Clair	DATE OF INSPECTION 05/27/2014
LOCATION OF INSTRUMENT (STREET AND CITY) 1 Paul Parks Dr. St. Clair		TIME OF INSPECTION 13:03 CDT

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or is operating within established limits. (Write in observed values where determined). Unmarked items must be corrected before using instrument.

☒ DIAGNOSTIC RECORD

<input checked="" type="checkbox"/> BLANK CHECK	<input checked="" type="checkbox"/> CO2 CHECK
<input checked="" type="checkbox"/> FC 1 TEMP	<input checked="" type="checkbox"/> FLOW CHECK
<input checked="" type="checkbox"/> SRC TEMP	<input checked="" type="checkbox"/> FCB CHECK
<input checked="" type="checkbox"/> DET TEMP	<input checked="" type="checkbox"/> CRC COMP CHECK
<input checked="" type="checkbox"/> BT TEMP	<input checked="" type="checkbox"/> CRC CAL CHECK
<input checked="" type="checkbox"/> STD 2 TEMP	<input checked="" type="checkbox"/> PRINT TEST
<input checked="" type="checkbox"/> ETH CHECK	

BREATH ANALYZER ACCURACY STANDARDS

<input checked="" type="checkbox"/> SIMULATOR SOLUTION	<input type="checkbox"/> COMPRESSED ETHANOL-GAS MIXTURE	
<input checked="" type="checkbox"/> STANDARD SUPPLIER Guth	LOT# 13280	EXP. DATE 10/16/2015
<input checked="" type="checkbox"/> SIMULATOR TEMP (34°C ±0.2°C) 34°C +/- .2°	SIMULATOR S/N SD2740	SIMULATOR EXP DATE 05/13/2015

☒ CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)

Run three tests using a standard solution. All three tests must be within +5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)

- ☒ 0.10% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE  
☐ 0.08% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE  
☐ 0.04% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 0.099 g/210L	TEST 2 0.099 g/210L	TEST 3 0.100 g/210L
---------------------	---------------------	---------------------

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:

REFUSALS	0	0-.04	0	.05-.09	0	.10-.14	0	.15-.19	0	OVER .19	0
----------	---	-------	---	---------	---	---------	---	---------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER

SIGNATURE 	PRINT FULL NAME SIKES, JOHN
TYPE II PERMIT NUMBER 240251	TELEPHONE NUMBER ( 636 ) 629-1313
EXPIRATION DATE 05/19/2016	

RETURN COMPLETED REPORT TO THE:

Breath Alcohol Program, Missouri Department of Health and Senior Services,  
Southeast District Office, 2875 James Blvd, Poplar Bluff, MO 63901



## GUTH LABORATORIES, INC.

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

### CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13280 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 18, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1217% (w/vol) ethyl alcohol. The expiration date for this lot number is October 16, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President  
GUTH LABORATORIES, INC.

**NIST Traceability:**

Testing was conducted using Cerilliant Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



STATE OF MISSOURI  
DEPARTMENT OF HEALTH AND SENIOR SERVICES  
BREATH ALCOHOL PROGRAM

2

**PERMIT**  
**TYPE II**  
**JOHN R SIKES**

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

**ALCO-SENSOR IV WITH PRINTER, INTOX EC/IR II**


for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.

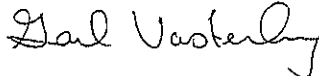
DATE 5/19/2014

NUMBER 240251

EXPIRES 5/19/2016

MO 580-0771 (8-10)

  
DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

  
DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

LAB-4 (R5-10)

	STATE OF MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES BREATH ALCOHOL PROGRAM
	<b>INSTRUMENT OPERATOR CARD</b>
<small>The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.</small>	
	
Operator SIKES, JOHN Permit No 240251 Date Issued 5/19/2014 Date Expires 5/19/2016	