



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES
STATE PUBLIC HEALTH LABORATORY
BREATH ALCOHOL PROGRAM
DATAMASTER MAINTENANCE REPORT

RECEIVED
By Carol Day at 8:11 am, Oct 08, 2014

REPORT #6

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).
Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.
Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

DATAMASTER SN 204081	NAME OF AGENCY MSHP	DATE OF INSPECTION 10/06/2014
LOCATION OF INSTRUMENT (STREET AND CITY) Scott County Sheriff's Office-211 S. New Madrid St.-Benton		TIME OF INSPECTION 12:36 pm

CHECKLIST: Place a mark in the box by each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	DATE AND TIME (from printout) 10/06/2014 12:38
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 50.0 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER
<input checked="" type="checkbox"/> INDICATOR LIGHTS	
<input checked="" type="checkbox"/> SIMULATOR SOLUTION SUPPLIER <u>Guth Laboratories Inc</u> LOT # <u>13290</u> EXP. DATE <u>10/29/2015</u>	
<input checked="" type="checkbox"/> SIMULATOR TEMP (34°C ± 0.2°C) <u>34.0</u> °C SIMULATOR SN <u>1099</u> EXP. DATE <u>02/18/2015</u>	
<input checked="" type="checkbox"/> CALIBRATION CHECK - (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)	
Run three tests using a standard solution. All three tests must be within ±5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)	
<input checked="" type="checkbox"/> 0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE	
<input type="checkbox"/> 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE	
<input type="checkbox"/> 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE	

TEST 1 • .098	TEST 2 • .100	TEST 3 • .100
---------------	---------------	---------------

☒ PERFORM R.F.I. TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:
(DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS	0	(0-.04)	0	(.05-.09)	0	(.10-.14)	0	(.15-.19)	1	OVER .19	0
----------	---	---------	---	-----------	---	-----------	---	-----------	---	----------	---

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS
(USE OTHER SIDE IF NECESSARY).

INSPECTING OFFICER

SIGNATURE <i>James Cooksey Jr.</i>	PRINT FULL NAME James C. Cooksey Jr.
TYPE / PERMIT NUMBER / EXPIRATION DATE 204086 03/11/2016	TELEPHONE NUMBER (573) 840-9500

RETURN COMPLETED REPORT TO THE:

Breath Alcohol Program, MO Department of Health and Senior Services, Southeast District Office
2875 James Blvd.
Poplar Bluff, MO 63901



GUTH LABORATORIES, INC.

680 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-584-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 13290 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 31, 2013, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1202% (w/vol) ethyl alcohol. The expiration date for this lot number is October 29, 2015 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN122211-02 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

BAC DataMaster Evidence Ticket

--- ORGANIZED CASE ---

(PL 105-105) (10/97) (10/97)

[illegible]

Operator Signature

James Cook



STATE OF MISSOURI
DEPARTMENT OF HEALTH AND SENIOR SERVICES
BREATH ALCOHOL PROGRAM

2

**PERMIT
TYPE II**

JAMES C COOKSEY JR.

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field service and repairs, and operate the following breath analyzer(s):

DATAMASTER, INTOX DMT


for the determination of the alcoholic content of blood from a sample of expired air. Permit issued under the provisions of sections 577.020 through 577.041, RSMo and 306.111 through 306.119 RSMo.


DATE 3/11/2014

NUMBER 240086


EXPIRES 3/11/2016

MO 580-0771 (6-10)

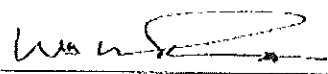
DIRECTOR OF STATE PUBLIC HEALTH LABORATORY

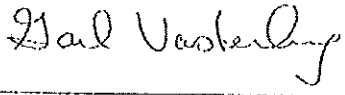

, acting director
DIRECTOR OF DEPARTMENT OF HEALTH AND SENIOR SERVICES

LAB-4 (116-10)

 STATE OF MISSOURI
DEPARTMENT OF HEALTH AND SENIOR SERVICES
BREATH ALCOHOL PROGRAM

INSTRUMENT OPERATOR CARD

The named cardholder is authorized to operate an evidential breath alcohol instrument for the determination of the alcoholic content in breath form of expired air in Missouri.



Operator COOKSEY JR., JAMES
Permit No 240086
Date Issued 3/11/2014 Date Expires 3/11/2016