



MISSOURI DEPARTMENT OF HEALTH AND SENIOR SERVICES  
STATE PUBLIC HEALTH LABORATORY  
BREATH ALCOHOL PROGRAM  
**CMI INTOXILYZER 5000 MAINTENANCE REPORT**

**RECEIVED**

By Carol Day at 10:04 am, Jul 14, 2014

Complete this report at the time of the regular monthly preventive maintenance check (not to exceed 35 days).  
Complete this report whenever the instrument is serviced or repaired and whenever it is placed into service.  
Retain the original and send a copy within 15 days to the Breath Alcohol Program, DHSS.

INTOXILYZER 5000 SN <u>66005190</u>	NAME OF AGENCY <u>Dexter Police Dept.</u>	DATE OF INSPECTION <u>7-3-14</u>
LOCATION OF INSTRUMENT (STREET AND CITY) <u>305 Cooper, Dexter MO 63841</u>		TIME OF INSPECTION <u>0700</u>

**CHECKLIST:** Place a mark by each item if found to be satisfactory or is operating within established limits. (Write in observed values where determined.) Unmarked items must be corrected before using instrument.

- ☒ DVM TEST: (.350 ± .150) .311
- ☒ DIAGNOSTIC CHECK (PRINTOUT ATTACHED) DATE AND TIME (FROM PRINTOUT) 7-3-14 0706
- ☒ CHARACTER DISPLAY TEST
- ☒ PRINT TEST (PRINTOUT ATTACHED)
- ☒ SIMULATOR SOLUTION SUPPLIER Repro LOT # 13002 EXP. DATE 6-9-15
- ☒ SIMULATOR TEMPERATURE (34°C ± 0.2°C) 34.0 C SIMULATOR SN SD3328 EXP. DATE 3-11-15
- ☒ CALIBRATION CHECK – (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)
- Run three tests using a standard solution. All three tests must be within ± 5% of the standard value and must have a spread of .005 or less. Mark the box corresponding to the standard solution being used. (PRINTOUT ATTACHED)
- ☐ 0.100% STANDARD - MUST READ BETWEEN 0.095% AND 0.105% INCLUSIVE
- ☐ 0.080% STANDARD - MUST READ BETWEEN 0.076% AND 0.084% INCLUSIVE
- ☐ 0.040% STANDARD - MUST READ BETWEEN 0.038% AND 0.042% INCLUSIVE

TEST 1 <u>.096</u>	TEST 2 <u>.096</u>	TEST 3 <u>.097</u>
-----------------------	-----------------------	-----------------------

☒ PERFORM RFI TEST (PRINTOUT ATTACHED)

INDICATE THE NUMBER OF BREATH TESTS IN THE FOLLOWING RANGES SINCE THE LAST MAINTENANCE REPORT:  
(DO NOT INCLUDE SELF-ADMINISTERED TESTS)

REFUSALS <u>0</u>	0-04 <u>0</u>	.05-.09 <u>0</u>	.10-.14 <u>2</u>	.15-.19 <u>0</u>	Over .19 <u>1</u>
-------------------	---------------	------------------	------------------	------------------	-------------------

LIST ANY NEW PARTS AND DESCRIBE ANY ALTERATION OR MODIFICATION THAT WAS MADE TO RESTORE THE INSTRUMENT TO OPERATE SATISFACTORILY AND WITHIN ESTABLISHED LIMITS (USE OTHER SIDE IF NECESSARY).

**INSPECTING OFFICER**

SIGNATURE <u>[Signature]</u>	PRINT FULL NAME <u>Sgt Joshua T. Benton</u>
TYPE II PERMIT NUMBER/EXPIRATION DATE <u>220243 9-11-14</u>	TELEPHONE NUMBER <u>573-624-5512</u>

RETURN COMPLETED REPORT TO THE: Breath Alcohol Program, Missouri Department of Health and Senior Services  
Southeast District Office  
2875 James Blvd.  
Poplar Bluff, MO 63901

**CERTIFICATE OF ANALYSIS**

**MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.**

**LOT NUMBER: 13002**

**EXPIRATION DATE: June 19, 2015 at 11:59 p.m.**

RepCo Marketing, Inc. certifies the following:


RepCo Marketing, Inc. manufactured, tested and supplied Lot Number 13002 of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed by an independent laboratory utilizing a gas chromatograph and found to contain .1217 gms/dl +/- .003 gms/dl wt/vol ethanol (95% Confidence).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of .100 +/- 3% gms/210L Breath when heated to 34 Degrees Celsius +/- 0.2 Degrees Celsius in a simulator (95% Confidence).


The date of manufacture for this lot number is June 20, 2013  
The expiration date for this lot number is June 19, 2015 at  
11:59 p.m.

This document is a true representation of the original Certificate of Analysis.

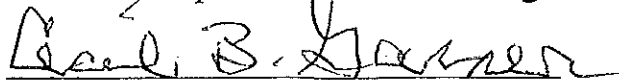

Cecil B. Garner, President  
RepCo Marketing, Inc.

State of Missouri  
DEPARTMENT OF HEALTH

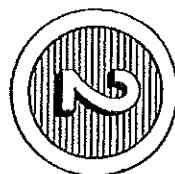

P E R M I T  
**TYPE II**

JOSHUA T BENTON

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field repairs, and operate the following breath analyzer(s):

INTOXILYZER 5000


for the determination of the alcoholic content of blood from a sample of expired (alveolar) air. Issued under the provisions of sections 577.020 through 577.041, RSMo 1986.


Date 09/11/2012

Number 220243

Expires 09/11/2014

MO 580-0771 (7-89)

  
Director of State Public Health Laboratory

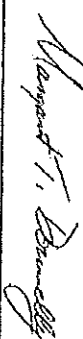  
Director, Department of Health

Lab. 4 (R7-88)

DEXTER PD  
INTOXILYZER - ALCOHOL ANALYZER  
NO MODEL 5000 SN 66-005190  
07/03/2014

DIAGNOSTIC TEST 07:06

PROM CHECK E735.23 PASSED  
RAM CHECK PASSED  
TEMP CHECK PASSED  
PROCESSOR CHECK PASSED  
SYNC PULSE PASSED  
SYNC SPEED PASSED  
NEG STABILITY PASSED  
POS STABILITY PASSED  
REF RANGE PASSED  
DIAGNOSTIC PASSED

PRINTER CHECK  
ABCDEFGHIJKLMNPQRSTUVWXYZ  
0123456789



INTOXILYZER® INSTRUMENT PRINTER CARD

0700  
Dexter PD  
Maintenance  
SUBJECT'S NAME  
INSTRUMENT LOCATION  
TIME FIRST OBSERVED  
OPERATOR  
ADDITIONAL INFORMATION AND/OR REMARKS

SN 66-005190  
E735.23  
INVALID TEST  
INHIBITED - RFI  
07/03/2014 07:06



INTOXILYZER® INSTRUMENT PRINTER CARD

0700  
Dexter PD  
Maintenance  
SUBJECT'S NAME  
INSTRUMENT LOCATION  
TIME FIRST OBSERVED  
OPERATOR  
ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER<sup>®</sup> INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

ADDITIONAL INFORMATION AND/OR REMARKS

OPERATOR

INSTRUMENT LOCATION

TIME FIRST OBSERVED

SUBJECT'S NAME

Maintenance



INTOXILYZER<sup>®</sup> INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

ADDITIONAL INFORMATION AND/OR REMARKS

OPERATOR

INSTRUMENT LOCATION

TIME FIRST OBSERVED

SUBJECT'S NAME

Maintenance

NO RFI PRESENT

TEST	%BAC	TIME
AIR BLANK	.000	07:14
CAL. CHECK	.096	07:15
AIR BLANK	.000	07:15
CAL. CHECK	.096	07:15
AIR BLANK	.000	07:15
CAL. CHECK	.096	07:15
AIR BLANK	.000	07:16

DEXTER PD  
INTOXILYZER - ALCOHOL ANALYZER  
NO MODEL 5000  
SN 66-005190  
07/03/2014

ABCEFGHIJKLMNOPQRSTUVWXYZ012345678910#1abcde  
ABCEFGHIJKLMNOPQRSTUVWXYZ012345678910#1abcde  
ABCEFGHIJKLMNOPQRSTUVWXYZ012345678910#1abcde  
ABCEFGHIJKLMNOPQRSTUVWXYZ012345678910#1abcde  
ABCEFGHIJKLMNOPQRSTUVWXYZ012345678910#1abcde  
ABCEFGHIJKLMNOPQRSTUVWXYZ012345678910#1abcde

SN 66-005190  
E735.23  
07/03/2014  
07:07